



CHANGING THE MIO® INFUSION SET: MINIMED® 630G INSULIN PUMP


START HERE:

1

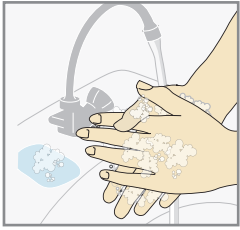


Wash your hands.
Press .

2

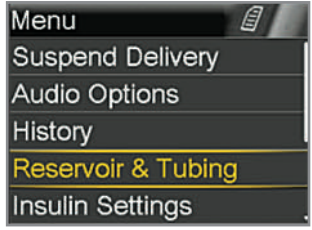


Select **Reservoir & Tubing**.

3





Select **New Reservoir**.

4

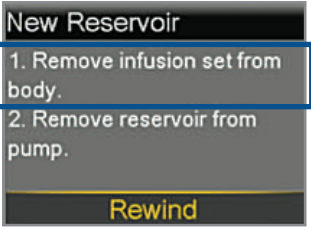


Remove the infusion set you have been using by loosening the adhesive and pulling away from body.

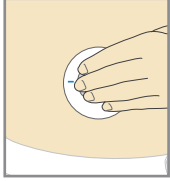


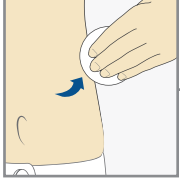


5

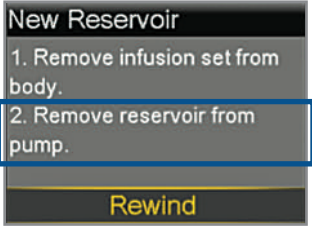


Remove the used reservoir from the pump.

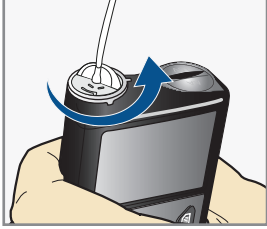

6




Select **Rewind**.




FILL RESERVOIR & CONNECT TO THE INFUSION SET TUBING
Follow the next steps to fill reservoir with insulin and connect to the infusion set tubing.

1

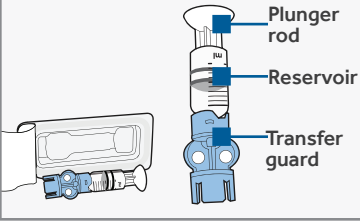


Remove from package.
Make sure insulin vial is at room temperature to reduce the risk of air bubbles.

2

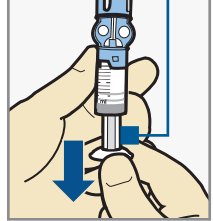


Pull plunger down to the amount that you plan to fill with insulin.

3

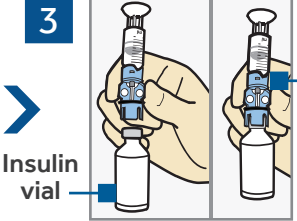

Wipe vial with alcohol. Place vial on table. Firmly press the blue transfer guard onto vial.

4

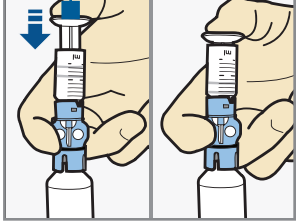

Push and hold plunger down.

Continued on next page

CHANGING THE MIO® INFUSION SET: MINIMED® 630G INSULIN PUMP

Warning: Do not use the reservoir or infusion set if any liquid gets on the top of the reservoir or inside the tubing connector. Liquid can temporarily block the vents and may result in the delivery of too little or too much insulin, which can cause hypoglycemia or hyperglycemia.

5

Insulin vial
Plunger

With your thumb still on the plunger, flip over so vial is on top. Release thumb and pull plunger down to fill with insulin.

6

Tap the reservoir to move air bubbles to top of reservoir. Push plunger up to move air into vial.

7

If needed, pull plunger back down to amount of insulin needed for 2-3 days.

8

Reservoir
Transfer guard

To avoid getting insulin on the top of the reservoir, turn vial over so it is upright. Hold transfer guard and turn reservoir counter-clockwise and remove from transfer guard.

CONNECT RESERVOIR TO INFUSION SET

You will place the tubing connector onto the end of the infusion set to the filled reservoir.

1

Remove the infusion set from packaging by pulling down the red plastic tab.

2

Pull the paper tab to remove the sterile paper backing from the bottom.

3

Squeeze the three raised markings on both sides of the lid. Lift and remove.

4

Free tubing from slot. Gently unwind tubing in a counter clockwise direction.

Connector

5

Gently push the tubing connector onto reservoir. Turn clockwise until locked. You will hear a click.

6

If you see air bubbles, tap reservoir to move them to top. Push plunger just a bit to move them into tubing.

7

Plunger

Twist plunger counter-clockwise to loosen and remove.

Continued on next page



CHANGING THE MIO® INFUSION SET: MINIMED® 630G INSULIN PUMP

THE BACKLIGHT MAY HAVE TURNED OFF

Press any button to turn the screen back on.




Select **Load Reservoir** and unlock pump if necessary.




Select **Next**.

PLACE RESERVOIR INTO PUMP


Now place the filled reservoir into the reservoir compartment of the pump.

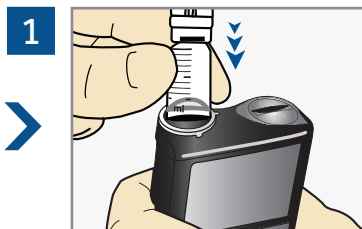

Place reservoir into pump.

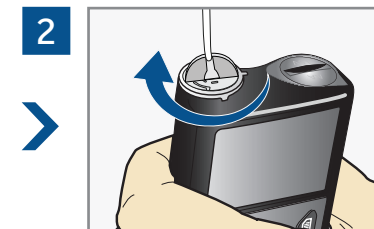


Turn clockwise until you feel reservoir lock into place.

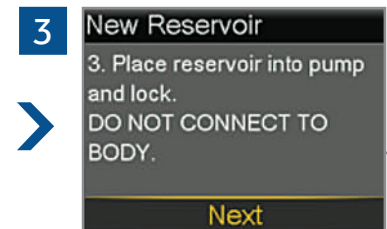


Select **Next**.

LOAD RESERVOIR AND FILL TUBING


Follow these steps to load the reservoir and fill the tubing.

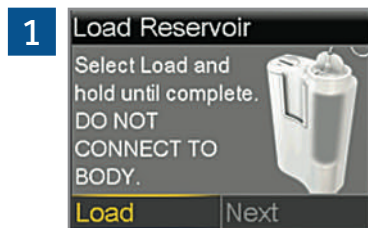


Select **Load** and keep holding .

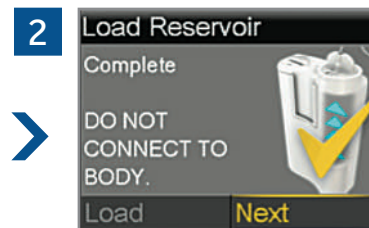

When you see this screen, select **Next**.

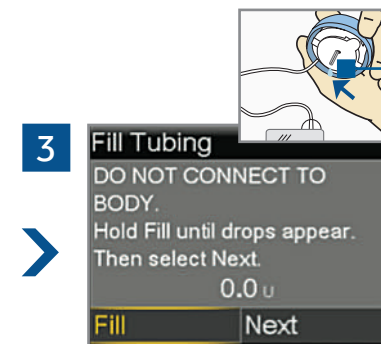

Select **Fill** and keep holding until you see drops at the end of tubing, then let go.

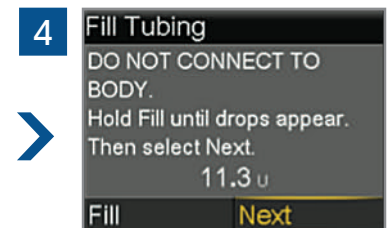

After you see drops, press and select **Next**.


Continued on next page




CHANGING THE MIO® INFUSION SET: MINIMED® 630G INSULIN PUMP

INSERT INFUSION SET

Next, follow the steps to insert the infusion set into your body.

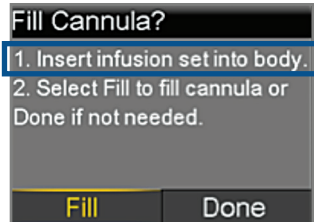


1

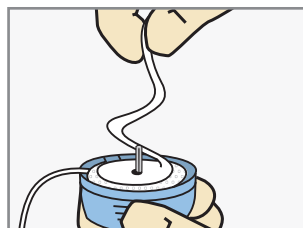


Gently peel paper to expose the adhesive.

2

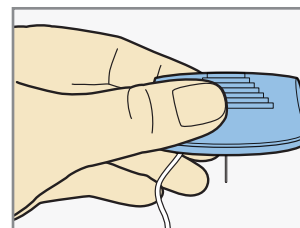


Turn over and hold by the lined ridges on the sides.

3

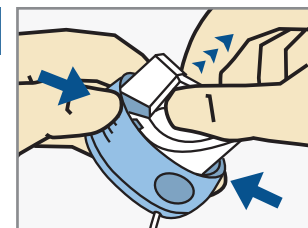


With other hand pull up on the center handle until it clicks and locks into place.

4

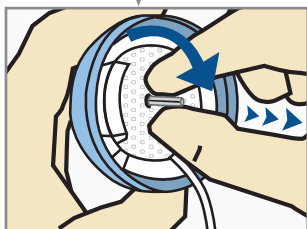

Carefully hold the needle guard by tip. Gently turn the needle guard and pull to remove it.

5

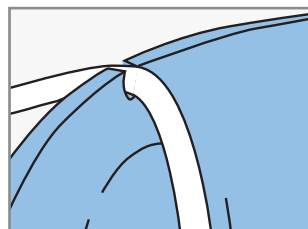


Place the tubing in slot on side of serter.



CAUTION:


Be careful not to pull the tubing tightly when placing it in the slot.

6

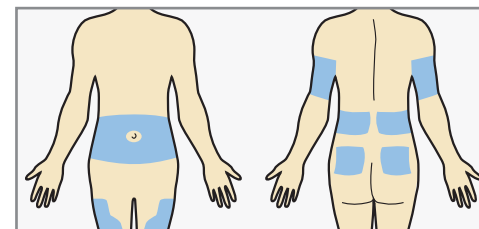


Choose an insertion site from the shaded areas shown here. Wipe with alcohol or other antiseptic.

7

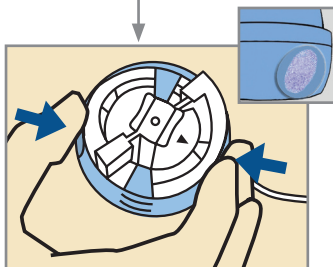


Place against prepared site on body. To insert, press the round indentations on each side of serter to insert needle.

8

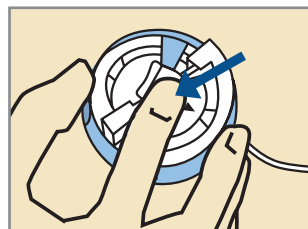

Push down on center of serter to press adhesive against skin.

9

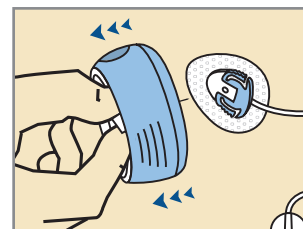

To remove serter, hold the center handle and gently pull it straight out away from body. Press the adhesive securely against skin.

10



Place the lid back onto serter to safely discard needle.


Continued on next page




CHANGING THE MIO® INFUSION SET: MINIMED® 630G INSULIN PUMP

FILL CANNULA

You will now fill the cannula, the little tube under your skin, with insulin.

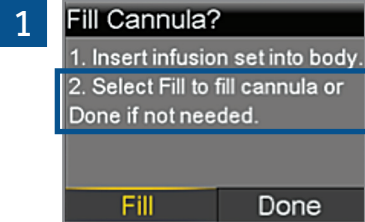

Select **Fill**.

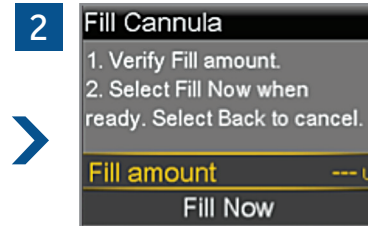

Select **Fill amount** and enter:

- 0.300 if using 6mm cannula
 - 0.500 if using 9mm cannula
- Then press **Fill Now**.




Note: Your pump will remember the **Fill amount** that you used last. Always verify that the Fill amount is correct.

- If it is correct, press **✓** to **Fill Now** and press **○**.
- If it is incorrect, press **○**. Change to correct amount and press **○**. Select **Fill Now**.



Select **Fill Now**.

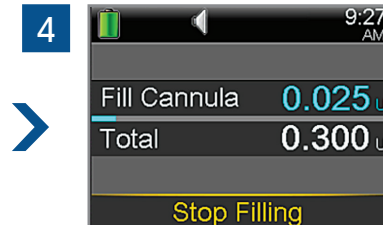

The Home screen displays the insulin as it fills the cannula.



Note: Select **Stop Filling** if you need to stop, for example, if you notice the **Total** amount is incorrect. This should rarely happen if you have verified the **Fill amount** on the previous screen.

**YOUR
INFUSION SET
CHANGE
IS NOW
COMPLETE!**